


Medtronic

Monitor longer. Fight secondary stroke.

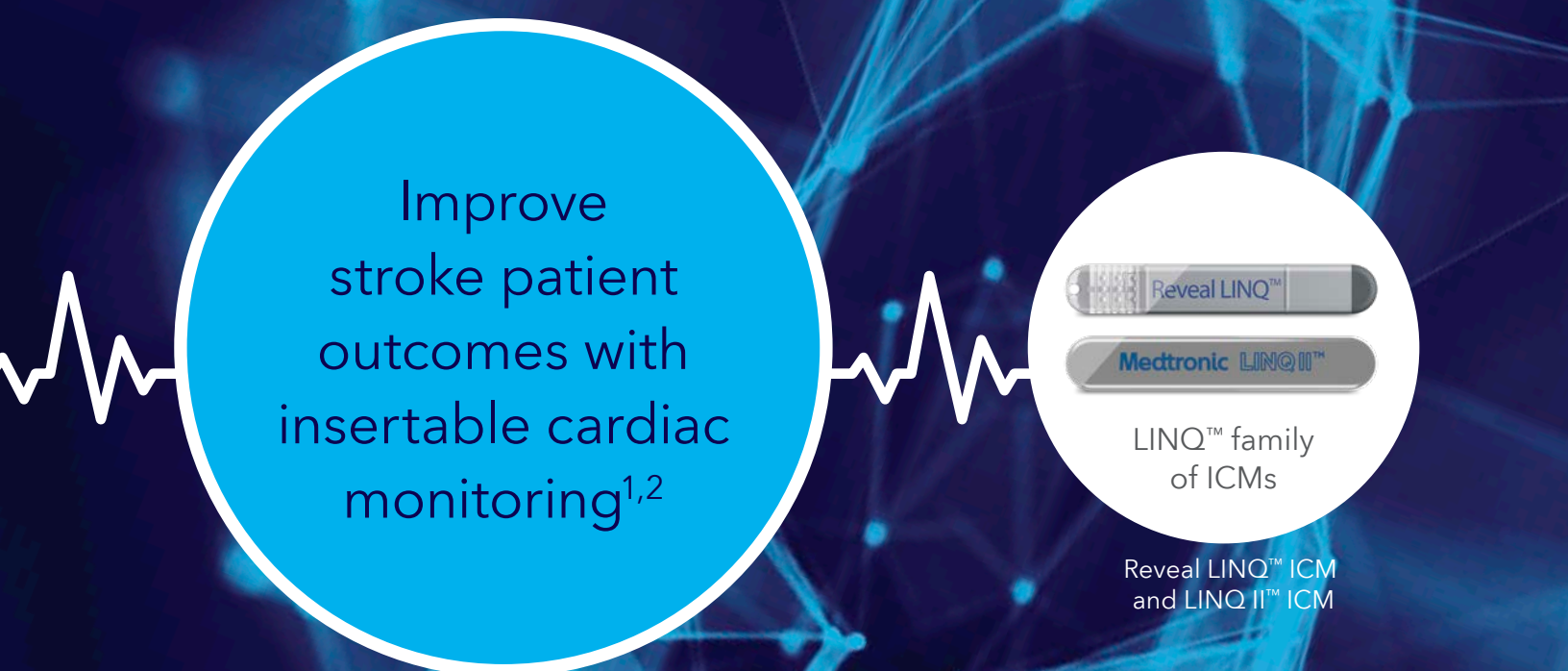

Improve
stroke patient
outcomes with
insertable cardiac
monitoring^{1,2}

Reveal LINQ™

Medtronic LINQ II™

LINQ™ family
of ICMs


Reveal LINQ™ ICM
and LINQ II™ ICM

Secondary stroke prevention is essential

Every year, more than

795,000

people in the United States have a stroke. Nearly 25% are in people who have had a previous stroke.³

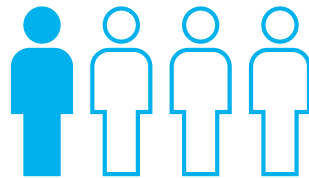

One in four stroke survivors will experience another stroke within five years.⁴




AF detection and treatment matter for improved patient outcomes.


Cryptogenic stroke and atrial fibrillation

Despite a comprehensive workup, 35% of ischemic strokes are cryptogenic.⁵

Cryptogenic stroke patients are at high risk for atrial fibrillation (AF). Up to 30% of patients with cryptogenic stroke have previously undetected AF up to three years post-stroke.⁶

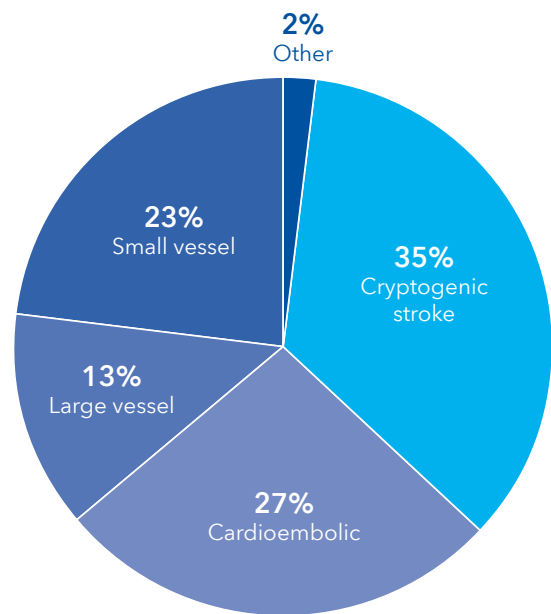

Detection of AF in cryptogenic stroke patients changes treatment⁵:

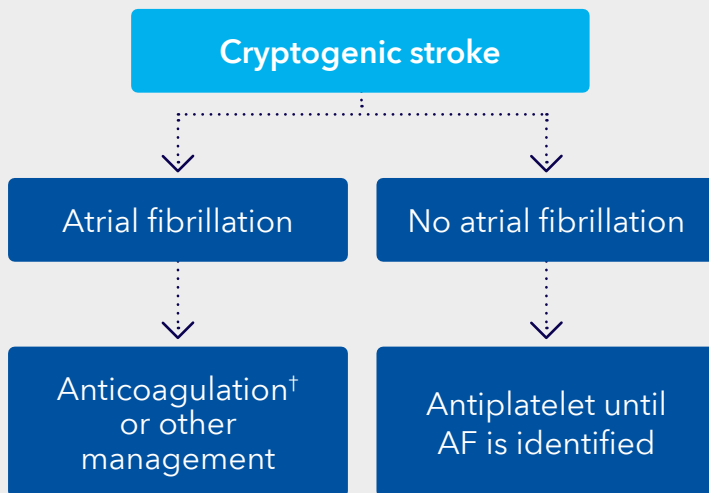


5x
AF is associated with a five-fold increase in the risk for ischemic stroke.⁷

2x
more likely for AF-related ischemic stroke to be fatal than non-AF stroke.⁸

†If patient is appropriate candidate.

Monitor longer

CRYSTAL-AF study results as published in the *New England Journal of Medicine*⁶

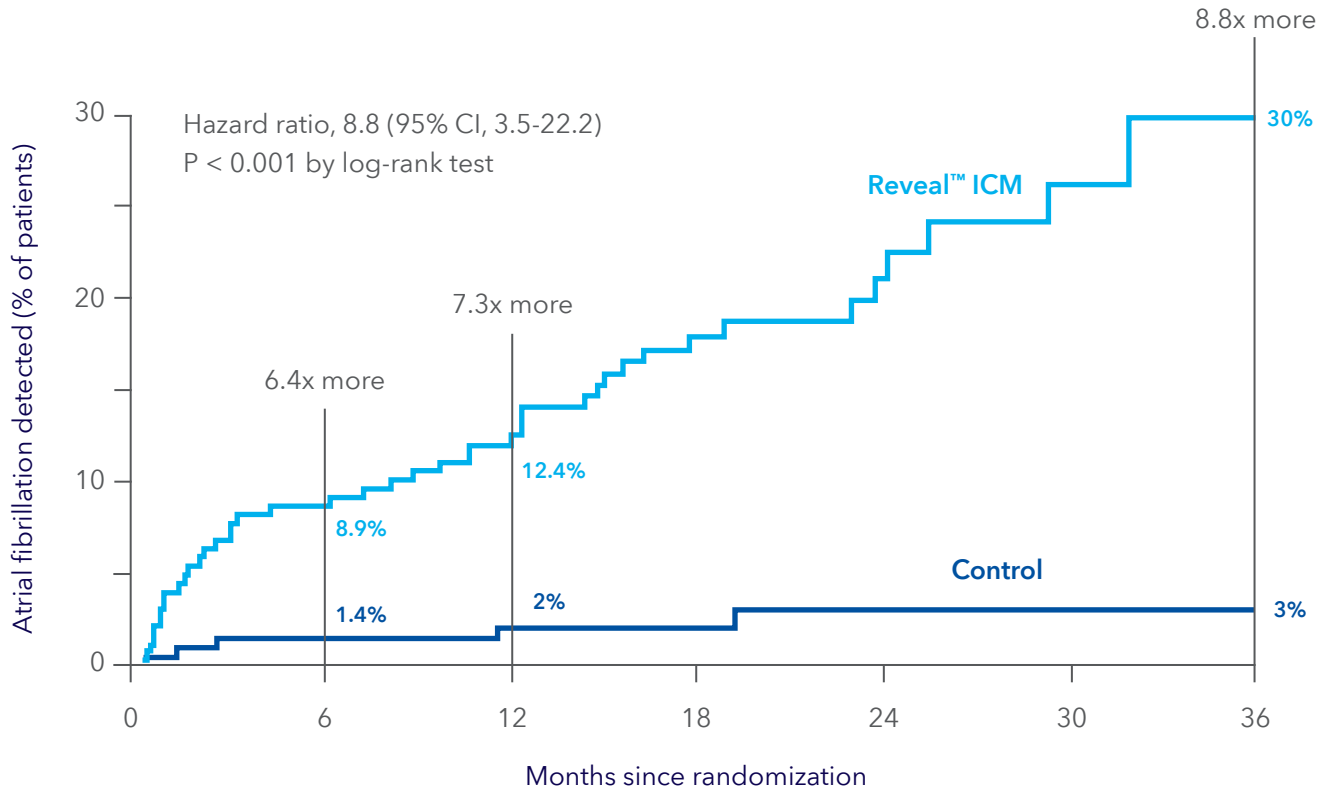

at risk

Control	220	194	167	114	72	36	7
ICM	221	191	173	102	57	29	8

88%

of patients who had AF would have been missed if only monitored for 30 days⁶

Based on Kaplan-Meier estimates.

79%

of first AF episodes were asymptomatic at 12 months⁶

97%


of patients with AF detected received oral anticoagulation⁶

Continuous monitoring with an insertable cardiac monitor (ICM) is superior to standard medical care for the detection of AF after cryptogenic stroke.⁶

External cardiac monitoring is not enough

Short-term monitoring could lead to undiagnosed patients⁶

Long-term, continuous monitoring (up to three years)

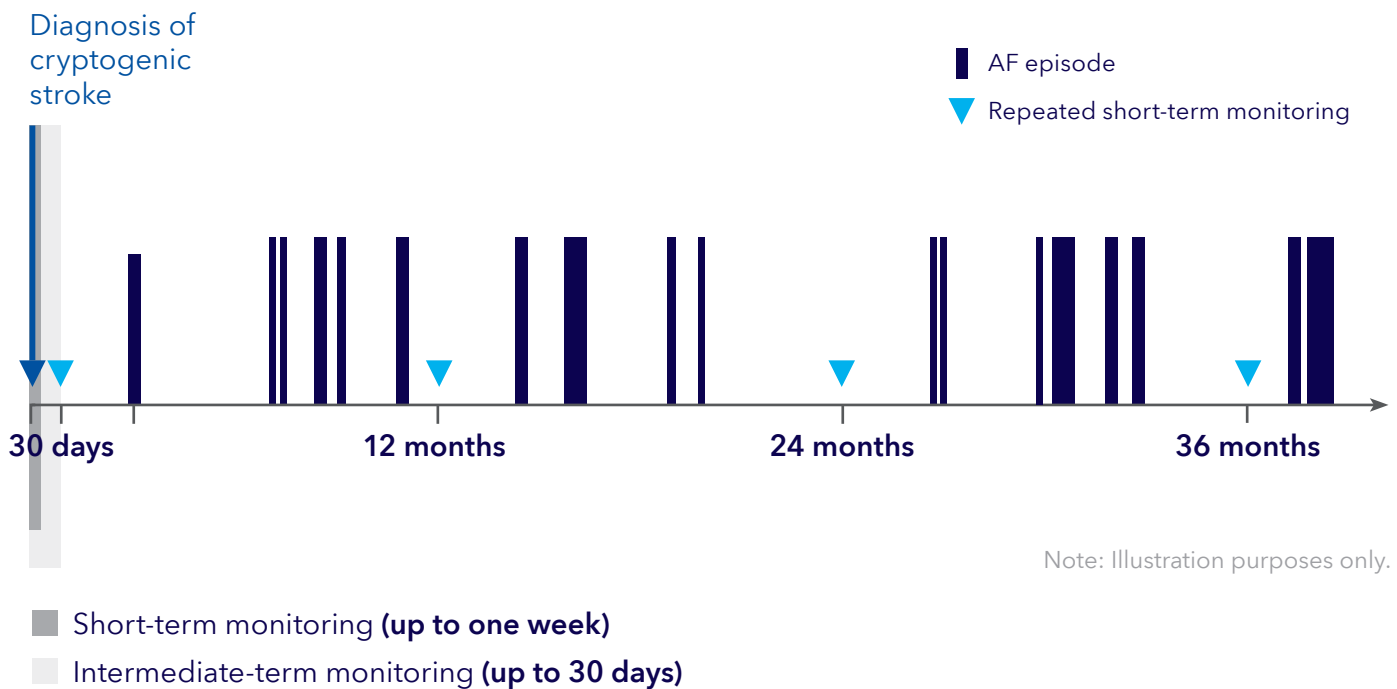

Medtronic ICM is superior to 30-day external loop recorder for AF detection⁹

15.3% AF detected in ICM arm ... versus ... **4.7%** in 30-day external loop recorder arm

Reveal LINQ ICM finds **3x** more AF over **12 months**.

Prevent secondary stroke

Cryptogenic stroke and TIA patients who underwent prolonged cardiac monitoring (PCM) compared to conventional cardiac monitoring show¹:

Find
AF.

Treat
AF.

Reduce
stroke.

2.5x

Increased
incidence of
AF detection.

2.1x

Increased
incidence of
anticoagulant
initiation.

55%

Decreased
risk of
recurrent
stroke.

Improve patient outcomes

Based on real-world claims data analysis^{2,10}



Faster time to AF diagnosis

Hazard ratio (95% CI) = 1.50 (1.40-1.60)



Faster time to OAC initiation

HR = 1.57 (1.42-1.73)



Reduced readmissions

21% decrease in annual readmissions
for ICM vs. ECM (p = 0.0006)



Reduced rate of mortality

HR = 0.70 (0.55-0.89)

Neurology and cardiology guidelines recommend ICM for cryptogenic stroke patients

2022

2022 ESO guideline for AF screening in adult patients with cryptogenic stroke and TIA¹¹

To maximize AF detection, clinicians should monitor cryptogenic stroke/TIA patients with implantable cardiac monitors, starting as soon as possible.

2021

2021 AHA/ASA guideline for the prevention of stroke in patients with stroke and TIA⁵

- Class 2a LOE B-R

In patients with cryptogenic stroke who do not have a contraindication to anticoagulation, long-term rhythm monitoring with mobile cardiac outpatient telemetry, implantable loop recorder, or other approach is reasonable to detect intermittent AF.

2020

2020 ESC guidelines for management of atrial fibrillation¹²

- Class 2a LOE B

In selected stroke patients without previously known AF, additional ECG monitoring using long-term, non-invasive ECG monitors or insertable cardiac monitors should be considered to detect AF.

2019

2019 AHA/ACC/HRS guideline for management of patients with AF¹³

- Class 2a LOE B-R

In patients with cryptogenic stroke (i.e., stroke of unknown cause) in whom external ambulatory monitoring is inconclusive, implantation of a cardiac monitor (loop recorder) is reasonable to optimize detection of silent AF.

Inform your clinical decisions with the LINQ family of ICMs

Reveal LINQ™ ICM and LINQ II™ ICM



Accuracy matters.

When you need answers, trust the most accurate insertable cardiac monitor (ICM) on the market.¹⁴⁻²⁶ Our proven ICMs deliver the data you need to make the right diagnosis for your patients, with continuous monitoring for up to 4.5 years.†

LINQ II ICM:

The world's most accurate ICM.¹⁴⁻²⁶

Reveal LINQ ICM:

The world's smallest ICM, powered by TruRhythm™ detection.²⁷

Industry-leading detection algorithms

The LINQ family of ICMs features the lowest published rate of AF false positives, and for greater confidence, the atrial fibrillation (AF) and Pause AccuRhythm™ AI algorithms further enhance the accuracy of the LINQ II ICM data.^{14-15, 28-30}

See what artificial intelligence can do.


The AccuRhythm™ AI platform is an artificial intelligence system that applies deep learning algorithms to LINQ II ICM data flowing into the CareLink™ network. The algorithms address the two most common sources of ICM false alerts – AF and pause.²⁸⁻³⁰

†Nominal settings

High patient satisfaction

More than twice as many patients reported being very satisfied with Reveal LINQ ICM versus the 30-day external loop recorder (43%, n = 54 in ICM arm, versus 20%, n = 26 in external loop recorder arm)⁹

- Medtronic ICMs are inserted just under the skin of the patient's chest in a short and simple procedure.
- The heart monitor is one-third the size of a AAA battery (1.2 cc) and is not visible in most patients.
- Use of Medtronic ICMs doesn't require a change in daily activities.

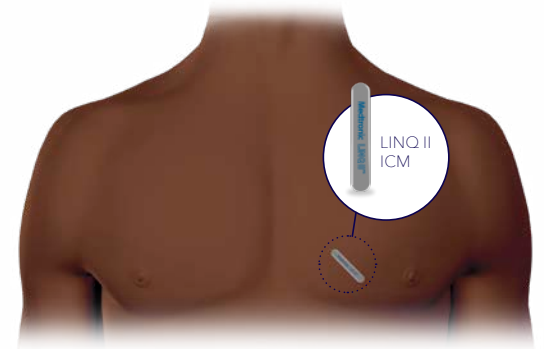

"It has put my mind, my heart, my body at rest knowing that I have the Medtronic LINQ."

- Pam
Reveal LINQ ICM patient

Scan QR code to watch Pam's story



Not every person will receive the same results.

Benefits of an ICM-first monitoring strategy are significant

- ✓ Greater AF detection yield with ICM^{6,9}
- ✓ Higher patient satisfaction with ICM⁹
- ✓ Patient compliance with short-term monitors is sub-optimal³¹
- ✓ Patients are lost to follow-up with short-term monitors^{2,32}
- ✓ Short-to-long monitoring causes delays in AF detection²
- ✓ ICM is cost-effective vs. SOC³³
- ✓ ICM is cost-saving vs. external monitors first³³


Less than 5% of ischemic stroke patients who initially receive a short-term monitor go on to receive an ICM.^{†2,32}

†Of patients with a clinical interaction with an electrophysiologist (EP) at any point in their care pathway.

Planning the cryptogenic stroke pathway



LINQ™ family of ICMs

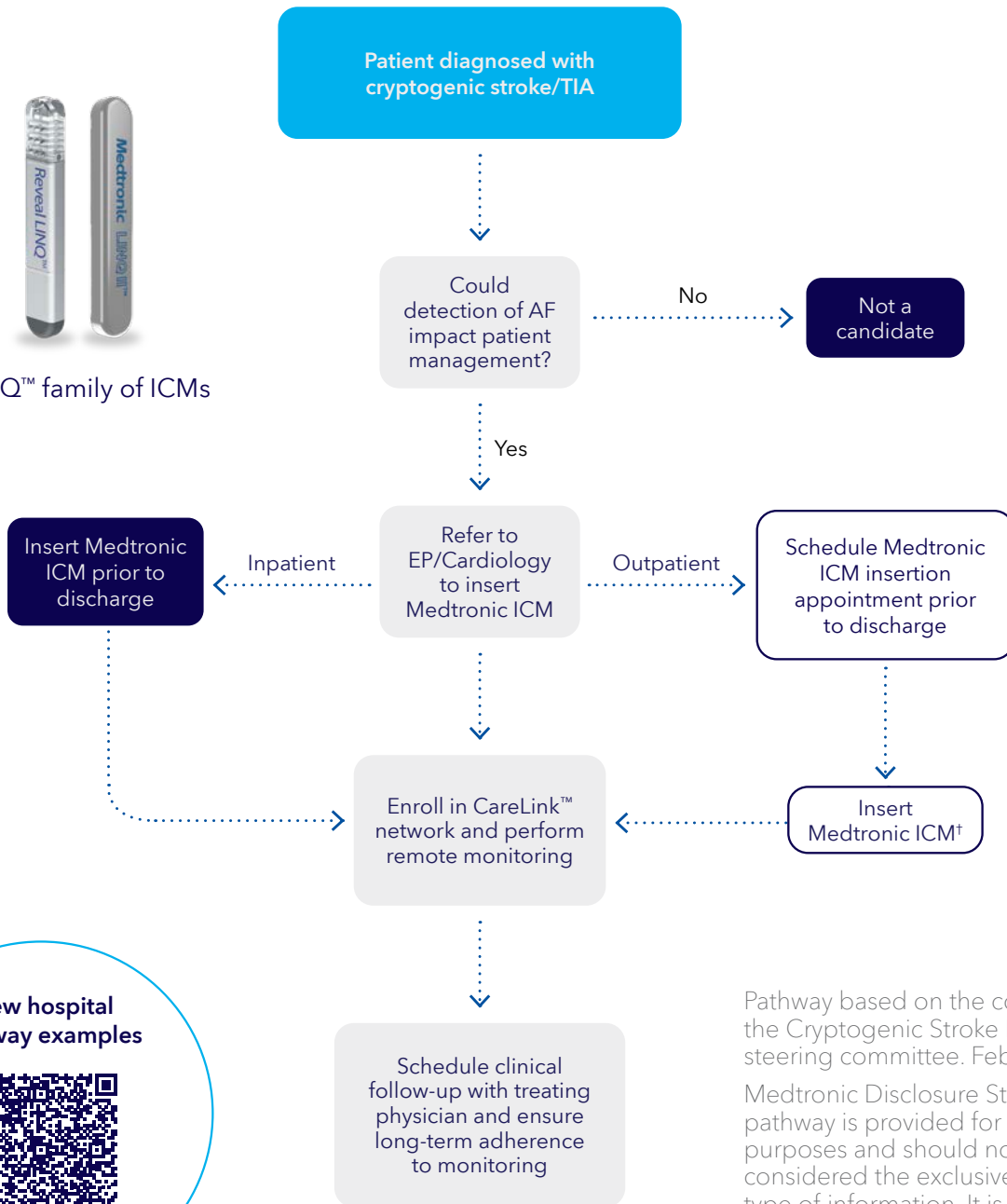

View hospital pathway examples



Pathway based on the consensus of the Cryptogenic Stroke Pathway steering committee. February 2016.

Medtronic Disclosure Statement: This pathway is provided for educational purposes and should not be considered the exclusive source for this type of information. It is the responsibility of the practitioner to exercise independent clinical judgment.

Refer to the brief statement for indications, warnings/precautions, and complications for Medtronic ICMs.

†If bridging with an external monitor before ICM, review short-term monitor results prior to ICM insertion. If AF is not detected, insert ICM.

Why do you need a cryptogenic stroke pathway?

Top reasons to pursue better multidisciplinary stroke care[‡]

For your patients



Better risk reduction strategy to prevent a secondary stroke



Ensures quality post-stroke care



Care focused on patient needs



Can help provide reassurance and peace of mind

For your hospital



Better risk reduction strategy to prevent a secondary stroke



Coordinated and integrated care



Consistent and reproducible approach



Enhanced hospital reputation for providing exemplary care



Stand out when pursuing Comprehensive Stroke Certification

[‡]Feedback based on the consensus of the Cryptogenic Stroke Pathway steering committee. February 2016.

ICM for large and small vessel stroke patients

STROKE AF study:


AF in non-cardioembolic stroke of presumed known origin.³⁴

Study design summary:

- Prospective, multisite, randomized clinical trial enrolling 496 patients at 33 centers in the United States
- Randomization 1:1 to continuous monitoring arm with Reveal LINQ ICM or control arm following site-specific standard of care (SOC) for detection of cardiac arrhythmias
- Follow-up: minimum 12 months, with maximum duration of 36 months


Primary objective

To determine whether long-term cardiac monitoring is superior to SOC for AF detection in patients with stroke attributed to large- or small- vessel disease through 12 months of follow-up.



45.5%
of first recurrent strokes are of a different etiology than the index stroke.³⁵

Detection of AF at 12 months

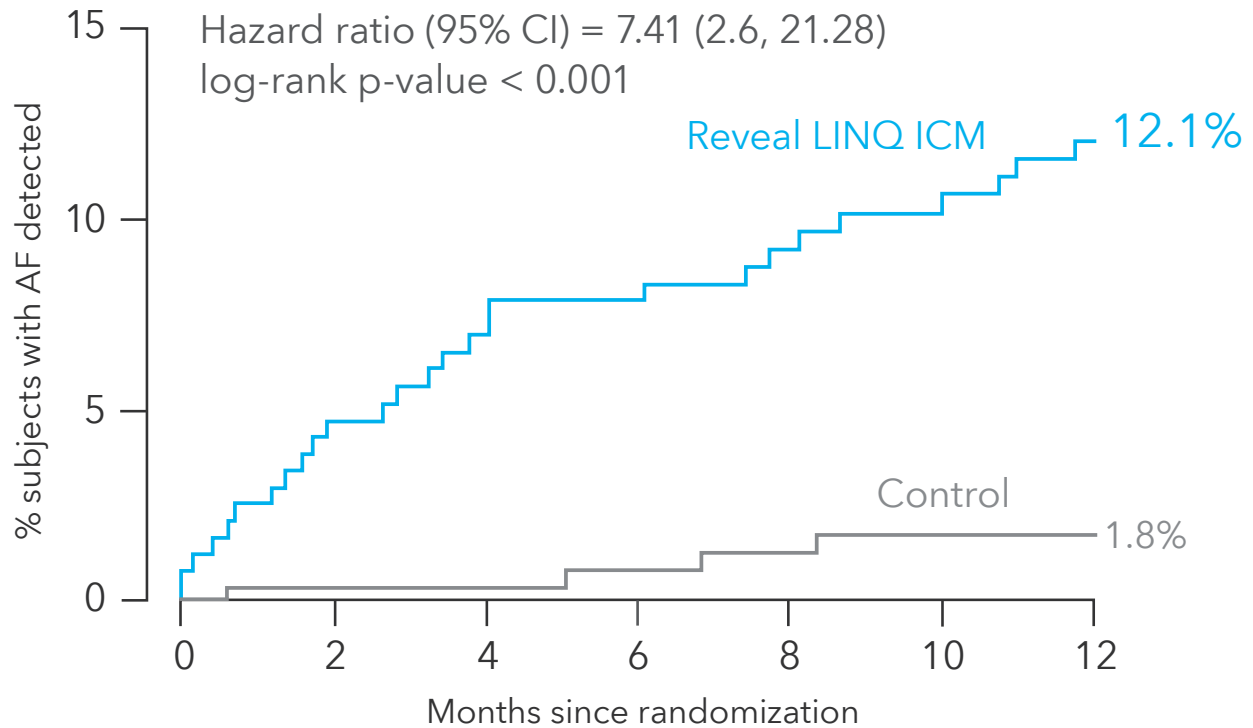

Results:

- Reveal LINQ ICM is superior to usual care for AF detection in large- and small-vessel stroke patients.
- 96.3% of first AF episodes were asymptomatic
- 55.5% of patients with AF detected with an ICM had an episode lasting greater than one hour.
- Median time to detection of AF was 99 days in the ICM arm.

Even the most aggressive intermittent monitoring strategy would have missed 78% of patients who had \geq two minutes of AF detected.

References

- 1 Tsivgoulis G, Katsanos AH, Grory BM, et al. Prolonged Cardiac Rhythm Monitoring and Secondary Stroke Prevention in Patients With Cryptogenic Cerebral Ischemia. *Stroke*. August 2019;50(8):2175-2180.
- 2 Yaghi S, Ryan MP, Gunnarsson CL, et al. Longitudinal outcomes in cryptogenic stroke patients with and without long-term cardiac monitoring for atrial fibrillation. *Heart Rhythm* O2. February 13, 2022;3(3):223-230.
- 3 Tsao CW, Aday AW, Almarazooq ZI, et al. Heart Disease and Stroke Statistics—2022 Update: A Report From the American Heart Association. *Circulation*. February 22, 2022;145(8):e153-e639.
- 4 Mohan KM, Wolfe CD, Rudd AG, Heuschmann PU, Kolominsky-Rabas PL, Grieve AP. Risk and cumulative risk of stroke recurrence: a systematic review and meta-analysis. *Stroke*. May 2011;42(5):1489-1494.
- 5 Kleindorfer DO, Towfighi A, Chaturvedi S, et al. 2021 Guideline for the Prevention of Stroke in Patients With Stroke and Transient Ischemic Attack: A Guideline From the American Heart Association/American Stroke Association. *Stroke*. July 2021;52(7):e364-e467.
- 6 Sanna T, Diener HC, Passman RS, et al. Cryptogenic stroke and underlying atrial fibrillation. *N Engl J Med*. June 26, 2014;370(26):2478-2486.
- 7 Wolf PA, Abbott RD, Kannel. Atrial fibrillation as an independent risk factor for stroke: the Framingham Study. *Stroke*. August 1991;22(8):983-988.
- 8 Lin HJ, Wolf PA, Kelly-Hayes M, et al. Stroke severity in atrial fibrillation. The Framingham Study. *Stroke*. October 1996;27(10):1760-1764.
- 9 Buck BH, Hill MD, Quinn FR, et al. Effects of Implantable vs Prolonged External Electrocardiographic Monitoring on Atrial Fibrillation Detection in Patients With Ischemic Stroke: The PER DIEM Randomized Clinical Trial. *JAMA*. June 1, 2021;325(21):2160-2168.
- 10 Yaghi S, et al. Healthcare Utilization in Cryptogenic Stroke Patients With vs. Without Long-term Cardiac Monitoring for Atrial Fibrillation: Evidence from Real World Data. Presented at World Stroke Congress 2021.
- 11 Rubiera M, Aires A, Antonenko, et al. European Stroke Organisation (ESO) guideline on screening for subclinical atrial fibrillation after stroke or transient ischaemic attack of undetermined origin. *Eur Stroke J*. Published online June 3, 2022.
- 12 Hindricks G, Potpara T, Dagres N, et al. 2020 ESC Guidelines for the diagnosis and management of atrial fibrillation developed in collaboration with the European Association for Cardio-Thoracic Surgery (EACTS): The Task Force for the diagnosis and management of atrial fibrillation of the European Society of Cardiology (ESC) Developed with the special contribution of the European Heart Rhythm Association (EHRA) of the ESC. *Eur Heart J*. 2021;42(5):373-498.
- 13 January CT, Wann LS, Alpert H, et al. 2019 AHA/ACC/HRS Focused Update of the 2014 AHA/ACC/HRS Guideline for the Management of Patients With Atrial Fibrillation: A Report of the American College of Cardiology/American Heart Association Task Force on Clinical Practice Guidelines and the Heart Rhythm Society. *J Am Coll Cardiol*. July 9, 2019;74(1):104-132.
- 14 Pürerfellner H, Sanders P, Sarkar S, et al. Adapting detection sensitivity based on evidence of irregular sinus arrhythmia to improve atrial fibrillation detection in insertable cardiac monitors. *Europace*. November 1, 2018;20(FI_3):f321-f328.
- 15 Nölker G, Mayer J, Boldt L, et al. Performance of an Implantable Cardiac Monitor to Detect Atrial Fibrillation: Results of the DETECT AF Study. *J Cardiovasc Electrophysiol*. December 2016;27(12):1403-1410.
- 16 Confirm Rx™ ICM K163407 FDA clearance letter. 2017.
- 17 Confirm Rx ICM K182981 FDA clearance letter. 2019.
- 18 Jot Dx™ ICM K212206 FDA clearance letter. 2021.
- 19 Monitoring devices Merlin PCS help manual for SJM Confirm, Confirm Rx ICM, Jot Dx manual. 2021.
- 20 BIOTRONIK BioMonitor™ 2 technical manual. 2017.
- 21 BIOTRONIK BIOMONITOR III technical manual. 2020.
- 22 BIOTRONIK BIOMONITOR III m technical manual. 2020.
- 23 BIOTRONIK BIOMONITOR III. K190548 FDA clearance. 2019.
- 24 BIOTRONIK BIOMONITOR III m. K201865 FDA clearance. 2020.
- 25 Lux-Dx™ ICM K212206 FDA clearance letter. 2020.
- 26 Lux-Dx ICM user manual. 2020.
- 27 ICM Size Comparison Guide. Medtronic data on file. 2021.
- 28 Cheng YJ, Ousdigian KT, Koehler J, Cho YK, Kloosterman M. Innovative Artificial Intelligence Application Reduces False Alerts while Maintaining Perfect True Pause Sensitivity for Insertable Cardiac Monitors. Presented at HRS 2021.
- 29 Radtke A, Ousdigian KT, Haddad TD, Koehler JL, Colombowala IK. Artificial Intelligence Enables Dramatic Reduction of False Atrial Fibrillation Alerts from Insertable Cardiac Monitors. Presented at HRS 2021.
- 30 AccuRhythm clinician manual supplements M015316C001 and M015314C001.
- 31 Kamel H, Navi BB, Eljovich L, et al. Pilot randomized trial of outpatient cardiac monitoring after cryptogenic stroke. *Stroke*. February 2013;44(2):528-530.
- 32 Landman SR and Sarkar S. Characterization of cardiac diagnostic care pathways by indication and physician specialty in a real-world dataset of 314,554 patients. Presented at ESC 2019.
- 33 Sawyer LM, Witte KK, Reynolds, MR, et al. Cost-effectiveness of an insertable cardiac monitor to detect atrial fibrillation in patients with cryptogenic stroke. *J Comp Eff Res*. February 2021;10(2):127-141.
- 34 Bernstein RA, Kamel H, Granger CB, et al. Effect of Long-term Continuous Cardiac Monitoring vs Usual Care on Detection of Atrial Fibrillation in Patients With Stroke Attributed to Large- or Small-Vessel Disease: The STROKE-AF Randomized Clinical Trial. *JAMA*. June 1, 2021;325(21):2169-2177.
- 35 Hillen T, Coshall C, Tilling K, et al. Cause of stroke recurrence is multifactorial: patterns, risk factors, and outcomes of stroke recurrence in the South London Stroke Register. *Stroke*. June 2003;34(6):11457-1463.

AccuRhythm AI ECG Classification System Brief Statement

Intended Use

The intended use of the system is to reduce false positive cardiac arrhythmia episodes.

Contraindications

There are no known contraindications for AccuRhythm AI Models ZA400, ZA410, or ZA420.

Precaution

The AccuRhythm AI algorithms ECG classification system may incorrectly adjudicate a true positive episode as an AI false episode, causing that episode to be suppressed in the remote monitoring system.

See the device manual for detailed information regarding the intended use, contraindications, warnings, precautions, and potential complications / adverse events. For further information, call Medtronic Technical Services at (800) 328-2518 and/or consult Medtronic's website at www.medtronic.com.

Caution: Federal law (USA) restricts these devices to sale by or on the order of a physician.

Brief Statement for Medtronic LINQ Family Insertable Cardiac Monitor (ICM) System and Remote Monitoring

Indications

The Reveal LINQ ICM is an insertable automatically-activated and patient-activated monitoring system that records subcutaneous ECG and is indicated in the following cases:

- patients with clinical syndromes or situations at increased risk of cardiac arrhythmias
- patients who experience transient symptoms such as dizziness, palpitation, syncope, and chest pain that may suggest a cardiac arrhythmia

This device has not been tested specifically for pediatric use.

The LINQ II ICM is an insertable automatically-activated and patient-activated monitoring system that records subcutaneous ECG and is indicated in adult patients, and in pediatric patients who are at least 2 years old, in the following cases:

- patients with clinical syndromes or situations at increased risk of cardiac arrhythmias
- patients who experience transient symptoms such as dizziness, palpitation, syncope, and chest pain that may suggest a cardiac arrhythmia

Contraindications

There are no known contraindications for the insertion of the LINQ Family ICMs or their accessories. However, the patient's particular medical condition may dictate whether or not a subcutaneous, chronically inserted device can be tolerated.

Warnings and Precautions

Patients with a LINQ Family ICM should avoid sources of diathermy, high sources of radiation, electrosurgical cautery, external defibrillation, lithotripsy, therapeutic ultrasound and radiofrequency ablation to avoid electrical reset of the device, and/or inappropriate sensing as described in the Medical procedure and EMI Warnings, Precautions and Guidance Manual. MRI scans should be performed only in a specified MR environment under specified conditions as described in the LINQ II or Reveal LINQ ICM MRI Technical Manual.

Wireless accessories available for use with a LINQ Family ICM may experience connectivity or performance issues. See product manuals for details and troubleshooting instructions.

Potential Adverse Events

Potential adverse events from the LINQ Family ICM include, but are not limited to, device rejection phenomena (including local tissue reaction), device migration, infection, and erosion through the skin.

There are no known adverse events associated with the use of any LINQ Family ICM wireless accessories.

See the device manuals for detailed information regarding the implant procedure, indications / intended use, contraindications, warnings, precautions, and potential complications/adverse events. For further information, please call Medtronic at (800) 328-2518 (Technical Services), (800) 551-5544 (Patient Services), and/or consult Medtronic's website at www.medtronic.com.

Caution: Federal law (USA) restricts prescription devices to sale by or on the order of a physician.

Medtronic

710 Medtronic Parkway
Minneapolis, MN 55432-5604
USA
Tel: (763) 514-4000
Toll-Free: (800) 328-2518

medtronic.com

©2023 Medtronic. All rights reserved. Medtronic, Medtronic logo, and Engineering the extraordinary are trademarks of Medtronic.™™™ Third-party brands are trademarks of their respective owners. All other brands are trademarks of a Medtronic company.

UC201405580g EN
03/2023